


# Attends® Slip Regular



Attends - The system with 10 absorption levels for bladder weakness

Attends Slip Regular are a range of all-in-one pads designed to manage heavy urinary and/or faecal incontinence.

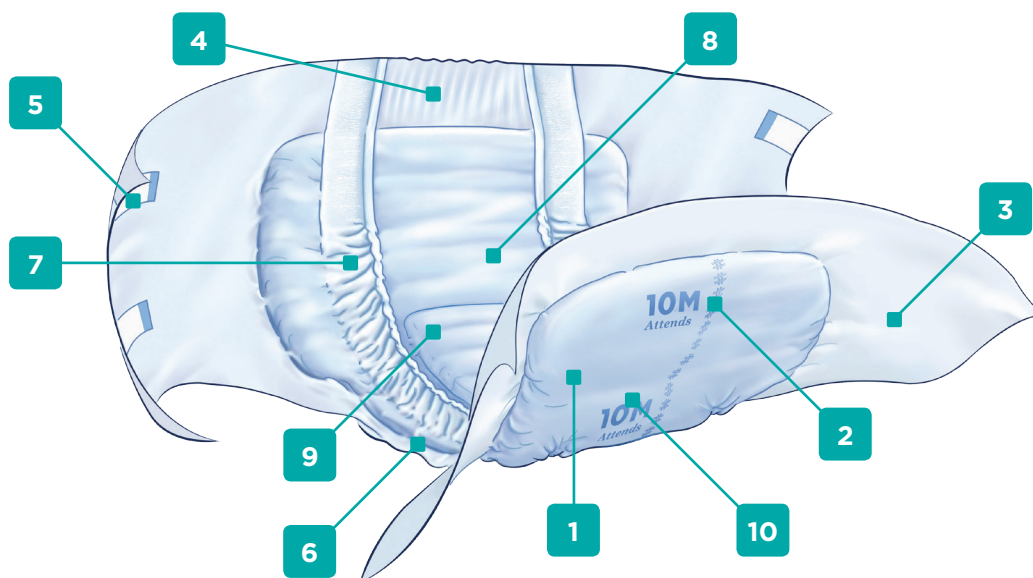

- |   |   |   |
|---|---|---|
| <p><b>1</b> White breathable textile like back sheet to aid containment except for size Extra Large.</p>                              | <p><b>2</b> Wetness indicator aids good nursing practice and acts as a guide to when the pad requires changing.</p> | <p><b>3</b> Breathable textile like wings provide comfort and support around the waist as well as improved skin health.</p>     |
| <p><b>4</b> Elasticated waistband for improved fit (Except for size XS).</p>  | <p><b>5</b> Re-sealable fixation system with four tapes to secure the pad in place.</p>                             | <p><b>6</b> Curved leg elastics for improved fit and to help reduce the risk of leakage.</p>                                    |
| <p><b>7</b> Full barrier leg cuff reduces the risk of leakage.</p>  | <p><b>8</b> Anatomically shaped core for a close comfortable fit.</p>   | <p><b>9</b> Super absorbent core with acquisition layer for improved leakage protection, skin dryness and odour protection.</p> |
| <p><b>10</b> The size and absorbency level are printed on the back sheet for easy identification of products in care homes/wards.</p> |   |   |

For more information visit [www.attends.com.au](http://www.attends.com.au)

	Packcount Casecount	EAN 7332152 Pack / Case	Hip / Waist Size (cm)	Backsheet	Absorption Core Material	Topsheet Material	Elasticated waistband	Leakage Barriers	Resealable fixation system	Breathable	Natural Odour Protection	Quick-Dry Technology
<b>Attends Slip Regular 8M</b>	28 4	203682 203699	80 - 110	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 8L</b>	28 4	203705 203712	100 - 150	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 9XS</b>	32 4	208182 208199	40 - 70	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✗	✓	✓	✓	✓	✓
<b>Attends Slip Regular 9S</b>	21 4	208106 208113	50 - 85	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 9M</b>	28 4	203316 203323	80 - 110	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 9L</b>	28 4	203330 203347	100 - 150	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 10S</b>	21 4	208120 208137	50 - 85	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 10M</b>	26 4	207680 207697	80 - 110	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 10L</b>	26 2	207659 207666	100 - 150	Laminate of nonwoven and breathable film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	✓	✓	✓
<b>Attends Slip Regular 10XL</b>	14 4	207635 207642	150 - 175	Laminate of polypropylene nonwoven and liquid impervious polyethylene film	Dual Core system with fluff and superabsorbents with Quick Dry acquisition layer	Non-woven	✓	✓	✓	Yes - Wings	✓	✓



**Quick-Dry**

Quick-Dry is the acquisition layer placed directly below the top sheet where the urine enters the pad. It quickly absorbs urine away from the skin and into the pad and prevents the urine returning to the surface. It can handle several gushes of urine without affecting how the pad works which results in exceptional skin dryness and leakage protection for the user



**pH Value**

The core technology with the curly fibre acquisition layer helps to guarantee a low pH which is similar to the skin. This helps to protect the skin and ensures improved skin health.



**Active Zone**

The Active Zone is an area centered in the middle of the pad, which provides fast absorption where it is really needed. This provides exceptional skin dryness and improved leakage protection for the user



**Natural Odour Protection**

Natural Odour Protection reduces the risk of unpleasant smells.



**Breathable**

The breathable material enables the skin to breathe, which contributes to improved skin health.



Attends Slips are tested and approved by proDERM, Institute for Applied Dermatological Research (does not include Slip Regular 9XS or 10XL)